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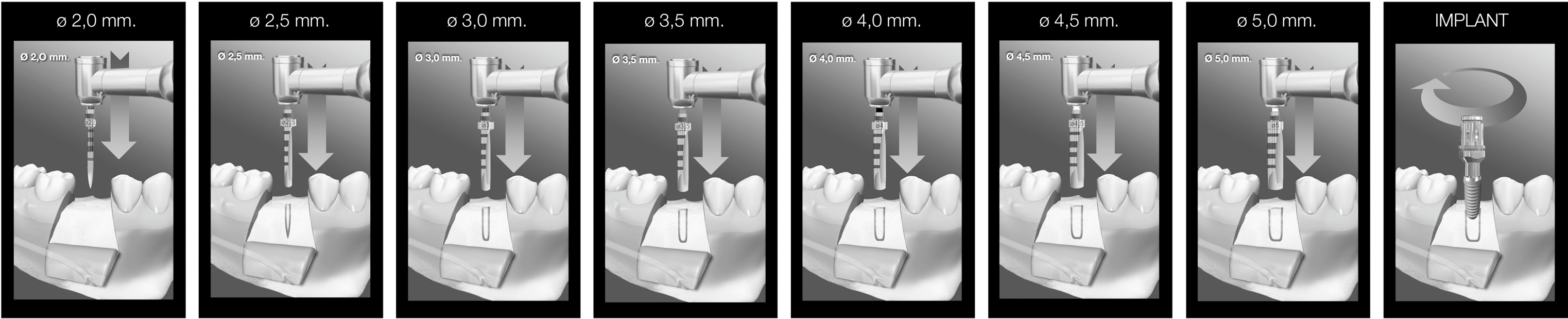
SURGICAL DRILLING SEQUENCE

After revealing the bone, start to drill with the Ø 2 lance drill, marking the site and direction. Drill down to the chosen intraosseous depth (equal to or longer than the length of the implant for placement). Revolutions: 1000-1200 rpm. Drill with the Ø 2.5 pilot drill to the chosen depth. The direction of the alveolus can be modified to adapt it to the planned implant axis. Revolutions: 1000-1200 rpm. Go through the various drill diameters until reaching the diameter of the chosen implant. Drill down to the height of the bed of the implant, at least, with each of the drills. Revolutions: 400 - 600 rpm. When drilling, check that irrigation is plentiful and that the cutting edge of the drills is in good condition. Drill with vertical movements in and out of the alveolus to cool the bone and prevent it from overheating. Eckermann's range of instruments includes drill stops that ensure that the drilling stops at the chosen depth, preventing further depth and involuntary contact with the anatomical structures of the jawbone. It is recommended the use of thread starters in situations of high density bone after milling, always at the discretion of the clinician. For bone regularization, Eckermann recommends the use of the ball mill (FBE 5). If the bone is of high density, after realization of the surgical bed the countersinking (FAE) can be used always under clinical criteria.

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